

bioretec

better healing - better life

ActivaNail™

- For fixation of bone fractures, osteotomies, arthrodesis and osteochondral fractures in the presence of appropriate immobilization.
- Self-Locking SL™ and grooved surface design enable nails to stay in place with improved rotational stability.
- ActivaNail™ is an implant with high-strength properties created by the material and the manufacturing technique.
- The surgeon has the option to slightly bend the implant in order to achieve the correct anatomical position.
- The bending modulus of the material is close to the value of cortical bone which decreases stress shielding in bone.
- Biodegrades safely within approximately 2 years with an excellent clinical track record.





Example of an indication: OCD (Osteochondritis Dissecans) with ActivaNail™



Manufacturer
Bioretec Ltd.
Hermiankatu 22, FI-33720 Tampere, Finland
Tel. +358 20 778 9500
sales@bioretec.com
www.bioretec.com

bioretec

Product codes and instruments

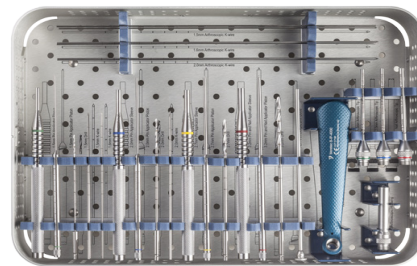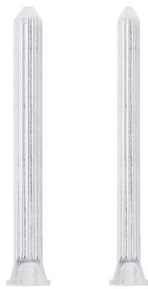

B-IPS-1000

- B-ANCO-1510 ActivaNail™ Conical 1.5 mm x 10 mm
- B-ANCO-1515 ActivaNail™ Conical 1.5 mm x 15 mm
- B-ANCO-1520 ActivaNail™ Conical 1.5 mm x 20 mm
- B-ANCO-1525 ActivaNail™ Conical 1.5 mm x 25 mm
- B-ANCO-1530 ActivaNail™ Conical 1.5 mm x 30 mm

- B-ANCO-2015 ActivaNail™ Conical 2.0 mm x 15 mm
- B-ANCO-2020 ActivaNail™ Conical 2.0 mm x 20 mm
- B-ANCO-2025 ActivaNail™ Conical 2.0 mm x 25 mm
- B-ANCO-2030 ActivaNail™ Conical 2.0 mm x 30 mm

B-IPS-1000 Instrument set, for ActivaPin™
(including instruments for ActivaNail™)

- B-IP-1500 Applicator Ø1.5 mm, reusable, for 1.5 mm pins
- B-IP-2000 Applicator Ø2.0 mm, reusable, for 2.0 mm pins
- B-IP-2700 Applicator Ø2.7 mm, reusable, for 2.7 mm pins
- B-IP-3200 Applicator Ø3.2 mm, reusable, for 3.2 mm pins
and for 1.5/2.0/2.7/3.2 mm nails
- B-IP-1501 K-Wire Ø1.6 mm, length 130 mm, for 1.5 mm pins/nails
- B-IP-1502 K-Wire Ø1.5 mm, length 130 mm, for 1.5 mm pins/nails
- B-IP-1506 K-Wire Ø1.5 mm, length 260 mm, for 1.5 mm pins/nails
- B-IP-2001 K-Wire Ø2.0 mm, length 130 mm, for 2.0 mm pins/nails
- B-IP-1504 K-Wire Ø1.6 mm, length 260 mm, arthroscopic,
or 1.5 mm pins/nails
- B-IP-2004 K-Wire Ø2.0 mm, length 260 mm, arthroscopic,
for 2.0 mm pins/nails
- B-IP-1503 Drill bit Ø1.5 mm, for 1.5 mm pins/nails, quick coupling
- B-IP-2003 Drill bit Ø2.0 mm, for 2.0 mm pins/nails, quick coupling
- B-IP-2703 Drill bit Ø2.7 mm, for 2.7 mm pins/nails, quick coupling
- B-IP-3203 Drill bit Ø3.2 mm, for 3.2 mm pins/nails, quick coupling
- B-IP-4000 Applicator handle, arthroscopic, for pins and nails
- B-IP-4001 Applicator drill guide, arthroscopic, for 1.5/2.0 mm nails
- B-IP-1512 Applicator tip, arthroscopic, for 1.5 mm pins
- B-IP-2012 Applicator tip, arthroscopic, for 2.0 mm pins
- B-IP-3212 Applicator tip, arthroscopic, for 1.5/2.0 mm pins
- B-IP-4004 Applicator housing, arthroscopic, for pins and nails
- B-IP-4003 Sterilisation tray, for ActivaPin™ instruments

REMARK: Prior to using Bioretec Implants examine thoroughly the Instructions For Use - inside each product package.

ORDERING INFORMATION: For more information please contact your nearest distributor or our customer service at Sales@bioretec.com.
Valid product range available from our website: www.bioretec.com.